

## The **DICOM virtual printer** for integrating all medical documents into the PACS



**Easily integrate your medical documents into the imaging network.**

**Print-in** makes it very simple to convert to DICOM format and to send all types of printable digital documents (PDF, Word, JPEG...) to the PACS.


Print-in software installs on any computer terminal (secretary, radiologist, physician, etc.) and acts like a **virtual printer**.

**Getting started** with Print-in is **very simple** and requires minimal user training.

Print-in is the **ideal and universal solution** for **integrating your medical documents** (radiology, ophthalmology, surgery, cardiology reports...).

Using Print-in you can link **imaging studies and supporting documents**, and consult them simultaneously.

Print-in is also the solution for integrating your reports on image CDs.

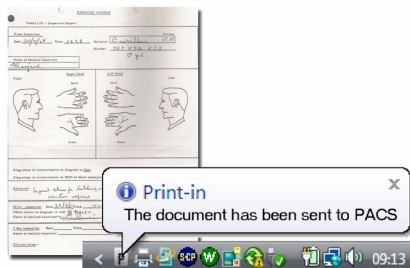

Minimal interface and maximum efficiency.

## Minimum System Requirements


Operating system: Windows 2000, XP, Vista

Any word processing or document digitizing software with Windows printing capability.



ETIAM Print-in

1. Select the "ETIAM - Print-in" printer.

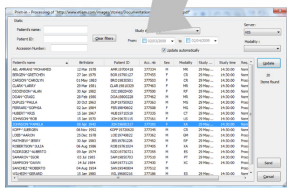

2. Select patient information from the HIS (DICOM Worklist) or from a study archived in the PACS.



3. Validate. Print-in converts the documents into DICOM and sends them to the PACS.

## Features

### Guaranteed compatibility

**Input:** All Windows software with printing capability.

**Output:** Conversion into DICOM PDF, DICOM SC (Secondary Capture, image or multiframe) for sending documents to the PACS or a DICOM printer.

### Easy start-up and storage management

Choice of DPI resolution for control of quality and size of files produced.

Simplified user interface, automation of each step. Requires no technical training.

### A plus for productivity

Integration with the information system: Use the worklist or a study archived in the PACS to assist patient selection.

Send documents on the imaging network as a background job.

## Companion Products

A natural complement to ETIAM products:

- Document, image, and video digitizing (DICOM Izer)
- DICOM report and image recording on CD/DVD (MARS)
- Viewing and selecting DICOM images and reports (Open LiteBox)
- Generating worklists in medical imaging (IDeal Broker)
- Storing, routing (STaR) or exchanging (SMN Router) images and reports

## Benefits

Getting started with Print-in is very simple and requires minimal user training.

Print-in improves the transfer of documents within the institution by making it possible to integrate all kinds of medical documents into the imaging network from a PC with Windows printing capability.

Print-in is compatible with all PACS on the market using a DICOM connection.

Use of the imaging network (PACS) infrastructure support improves access to documents and reduces printing costs.

Print-in enables to link imaging studies and supporting documents, and consult them simultaneously.

Solution for integrating reports on image CDs (MARS solutions).

Print-in transfers documents associated with other imaging studies intended for other institutions by means of a teleimaging system (SMN Router solutions).

This document is not contractual and this offer is subject to changes without prior notice.

All trademarks and registered trademarks herein are the property of their respective owners.

Distributed by: